



Misty Max 10™  
Disposable nebulizer

## Misty Max 10™ disposable nebulizer

An advanced jet design for:

- Fast treatment times
- High percentage of respirable particles
- Low residual volume
- Consistent and reliable performance

### Top 10 Features

- 1** Full 10 cc capacity provides medication dosing flexibility.



- 2** Unique, integrated (jet/baffle) design increases consistency and reliability in medication delivery.

- 3** Tinted snap-in jet remains in place yet allows for removal during cleaning.

- 4** Single twist closure of top and bottom gives confirmation of proper assembly.

- 5** Secure top and bottom seals prevent leakage.

- 6** Anti-spill design helps prevent medication waste.

- 7** Angled operation accommodates various patient positions.

- 8** Jet technology delivers efficient, high output rates, resulting in fast treatment times.

- 9** Clear base with visible graduations facilitates medication loading and visualisation of aerosol delivery.

- 10** Easy to assemble rigid tubing connector and crush-resistant tubing.






Delivering fast and effective treatment, more efficiently. Accommodating all types of patients in all types of situations. These were the forces behind the development of the Misty Max 10™ nebulizer. Featuring a new type of jet design, the Misty Max 10™ nebulizer can speed up treatment and help ensure consistent medication delivery while offering you more flexibility in treatment options.

Nebulizer is available preassembled and fits virtually all common masks and adapters. All oxygen tubing is crush-resistant. Recommended for single patient use. The one-piece jet design and single-thread top connection make the unit easy to disassemble and reassemble.

**Technical information (@ 8 L/min flow)\***

MMAD	1.61 μm
GSD	2.18
Respirable fraction % < 4.7μm (by mass)	85.2
Output rate	0.325 g/min
Treatment time	8 min, 9 sec
Residual volume	0.92 g
Fill capacity	10 mL

\* Please see the Misty Max 10 nebulizer technical data sheet for additional information, including performance characteristics at 5 L/min and 10 L/min flow.

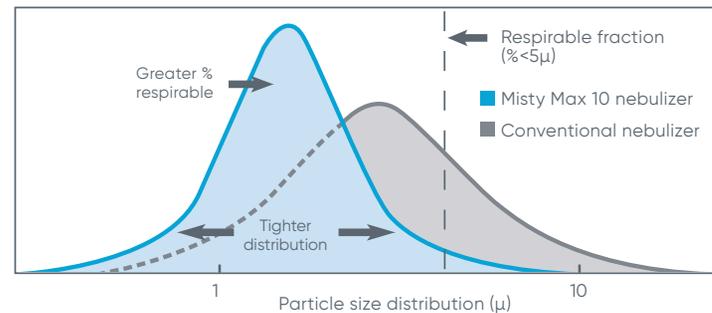

Particle size is an important characteristic of nebulizer performance. The Misty Max 10 nebulizer provides fine particles with a tighter distribution.\*

\* As proven by in-house testing

## Ordering information

Cat. no.	Description
Nebulizers without mask, 50 per package	
002430	Nebulizer only
002431	With 7' oxygen tubing
002434	With baffled tee adapter, 7' oxygen tubing and mouthpiece
002435	With tee adapter, 7' oxygen tubing and mouthpiece
002438	With baffled tee adapter, 7' oxygen tubing, mouthpiece and 6" flextube
002446	With baffled tee, 7' U/Connect-It oxygen tubing, mouthpiece and 6" flextube
Nebulizer without mask, 25 per package	
002450	With patient "Y" mouthpiece, 7' oxygen tubing, 6" flextube and bacteria filter

Cat. no.	Description
Nebulizers with adult mask, 50 per package	
002433	With adult aerosol mask and 7' oxygen tubing
Nebulizers with pediatric mask, 50 per package	
002443	With pediatric aerosol mask
002444	With pediatric aerosol mask and 7' oxygen tubing
Masks only, 50 per case	
001261	Pediatric aerosol mask, short style
001263	Pediatric aerosol mask, under-the-chin style
001206	Adult aerosol mask

For more information, please contact your local Vyaire representative.



Telephone: 02-86359747, 02-86379464

Fax: 02-86359748

e-mail: [sales@lifelink-inc.com](mailto:sales@lifelink-inc.com)